

“Designed
for Superior
Results.”

simple · powerful · flexible

The highly efficient fine-blasting devices of the Cemat series guarantee top-quality results. These robust machines, with stainless steel housing and 2 to 4 full-metal blasting medium containers, offer everything you would expect from a device of this class – high comfort included!

OUR PROVEN FORMULA FOR SUCCESS:

We combine technology, experience, expertise, and team spirit.



Always here for you

Do you want to request a quote, order products, or have questions about your order, delivery, or warranty?

We look forward to your call and are happy to assist you.



www.wassermann.hamburg



fb.com/WassermannDentalmaschinen

**Cemat &
Cemat Plus**



Fig.: Cemat-4

Info Hotline

+49 (0)40 730 926-0



SCAN ME


Cemat-2 Plus Cemat-3 Plus Cemat-4 Plus

Cemat Plus:
With intuitive handle-trigger

QUALITY MADE IN GERMANY SINCE 1927

Whether for blasting, roughening, or polishing:

Details like the color-coded system and effective lighting of the work area with a protected LED light make the job easier. The entire system is always under pressure and thus immediately ready for use.

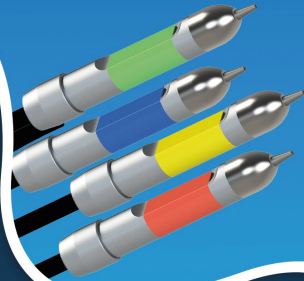

Key features at a glance:

- Versatile fine-blasting devices with 2 to 4 full-metal chambers
- Handy air-blow nozzle for easy cleaning
- Ultra-fast and uniquely precise start/stop technology
- Adjustable pressure and blast volume per chamber
- Saves costs with minimal consumption
- High hygiene standards with replaceable grip sleeves
- Specially durable blasting hoses
- Convenient sand drawer
 - Cemat-4 with 4 fine-blasting nozzles: 0.6 (yellow), 0.8 (blue), 1.2 (red), 0.8 (green)
 - Cemat-3 with 3 fine-blasting nozzles: 0.8 (blue), 1.2 (red), 0.8 (green)
 - Cemat-2 with 2 fine-blasting nozzles: 0.8 (blue), 1.2 (red)

Technical Specifications

Working Pressure: 0.8–6 bar
Power Consumption: 0.2 A

Dimensions W x H x D: 445 x 460 x 620 mm
Weight C-4, C-3, C-2: 28 kg, 26 kg, 24 kg

Sound Level: ≤ 70 dB(A)
Tank Capacity: 560 cm³

Info & Consultation:

Wassermann Dental-Maschinen GmbH
Rudorffweg 15 - 17
21031 Hamburg, Germany

Tel: +49 (0)40 730 926-0
www.wassermann-dental.com
sales@wassermann-dental.com